



**AT300 / AT600 Series**  
**Precise Body Temperature Measurement Thermal Camera**  
**Quick Start Guide V1.0.3**

IRay Technology Co., Ltd.

©IRay Technology Co., Ltd. 2018 Reserve all the right. All in this manual including texts, pictures, diagrams and other contents belong to IRay Technology Co.,Ltd (Hereinafter referred to as “Our company”or “IRay Photoelectric”). Without the written permission, no one shall copy, photocopy, translate or disseminate all or part of this manual.

This manual is used as a guide. The photos, graphics, diagrams and illustrations provided in the manual are only used to explain, which may be different from the specific product. Please refer to the real object. We try our best to make sure the contents in this manual are accurate. We do not provide any representations or warranties in this manual.

If you need the latest version of this manual, please contact us. IRay Photoelectric recommends that you use this manual under the guidance of professionals.

## Version History

| Version | Date       | Description                     |
|---------|------------|---------------------------------|
| V1.0.0  | 2020-03-06 | Initial version                 |
| V1.0.1  | 2020-03-11 | Modify the expression in page 9 |
| V1.0.2  | 2020-03-19 | Modify standard packing list    |
| V1.0.3  | 2020-03-25 | Modify the product name         |

## Contents

|  |        |
|--|--------|
| 1. Overview .....                                    | - 3 -  |
| 2. Product Performance .....                         | - 4 -  |
| 3. Installation Guide .....                          | - 6 -  |
| 3.1 Product and Accessories List .....               | - 6 -  |
| 3.2 Mounting and Connecting Thermal Imager.....      | - 7 -  |
| 3.3 Introduction of Reference Blackbody .....        | - 10 - |
| 3.3.1 Specification .....                            | - 10 - |
| 3.3.2 Structure of Reference Blackbody .....         | - 10 - |
| 3.3.3 Quick Start and Attentions .....               | - 11 - |
| 3.4 Site Layout Solution .....                       | - 12 - |
| 3.4 Site Installation and Layout Considerations..... | - 13 - |
| 4 Quick Start of Client .....                        | - 14 - |
| 5、 Troubleshooting .....                             | - 15 - |
| 6、 Maintenance .....                                 | - 16 - |
| 6.1 Notice .....                                     | - 16 - |
| 6.2 Cleaning and Maintenance .....                   | - 16 - |
| 6.3 Safety Codes .....                               | - 17 - |
| 7. After sales service .....                         | - 17 - |
| 8. Company information .....                         | - 17 - |

## 1. Overview

AT300 / 600 accurate temperature measurement thermal camera can be widely used in medical temperature measurement, body temperature screening, animal husbandry, machine vision, high-precision power detection, high-precision industrial detection and other industries. It can realize real-time temperature information transmission, accurate detection, fast screening and so on. In particular, its abnormal body temperature screening function can be widely used in airports, railways, stations, terminals, schools, and other densely populated places, to achieve rapid screening and avoid congestion.



Figure 1 Product appearance

## 2. Product Performance

Table 1 Product performance parameter









| Model   | AT3003F   | AT6005F     |
|---|---|-------------|
| <b>Performance</b>  |   |             |
| <b>Detector type</b>  | VOx Uncooled Infrared FPA Detector  |             |
| <b>Resolution</b>   | 384×288   | 640×512     |
| <b>Lens Model</b>   |   |             |
| <b>Focal Distance</b>   | 7.8mm   | 15mm        |
| <b>Lens control</b>   | Support auto focus  |             |
| <b>Focus Type</b>   | Athermal  | Athermal    |
| <b>FOV(Horizon×Vertical)</b>  | 47°×35.6°   | 34.2°×27.4° |
| <b>IFOV</b>   | 2.1mrad   | 0.93mrad    |
| <b>Temperature Measurement Function</b>   |   |             |
| <b>Measuring Range</b>  | 0℃~50℃  |             |
| <b>Measurement accuracy@ ambient temperature 16℃ ~ 32℃ (with blackbody)</b>       | ±2℃ @Target temperature 0℃~20℃<br>±0.6℃ @Target temperature 20℃~33℃<br>±0.3℃ @Target temperature 33℃~42℃<br>±0.6℃ @Target temperature 42℃~50℃   |             |
| <b>Effective temperature measurement distance (meet the measurement accuracy)</b> | 1m~5m   | 1m~5m       |
| <b>Temperature measurement tool</b>   | Arbitrary fixed spots<br>Full screen highest / lowest temperature spot capturing<br>Center spot temperature measurement<br>Line / area analysis tools<br>Manual temperature width selection |             |
| <b>Ethernet</b>   |   |             |
| <b>Network Interface</b>  | RJ45  |             |

|   |   |       |
|---|---|-------|
| <b>Network Protocol</b>                   | TCP、UDP、ICMP、IGMP、DHCP、RTSP                             |       |
| <b>Image Adjustment</b>                   |   |       |
| <b>Brightness and contrast adjustment</b> | Manual / Auto 0 (default) / Auto 1                      |       |
| <b>Polarity</b>                           | Black hot/White hot                                     |       |
| <b>Palette</b>                            | 18 modes are supportable                                |       |
| <b>Image flip</b>                         | Left -right / up-down / diagonal                        |       |
| <b>Power Supply</b>                       |   |       |
| <b>Supply voltage</b>                     | 10-36VDC  |       |
| <b>Typical consumption @25°C</b>          | ≤3W   | ≤3.3W |
| <b>Power protection</b>                   | Support   |       |
| <b>POE power supply</b>                   | Support (based on IEEE 802.3af)                         |       |
| <b>Physical Property</b>                  |   |       |
| <b>Size</b>                               | 55mm×55mm × 119mm (L × W × H)                           |       |
| <b>Environmental Adaptation</b>           |   |       |
| <b>Operating temperature</b>              | -10°C ~ 60°C (16~32°C accurate temperature measurement) |       |
| <b>Storage temperature</b>                | -20°C ~ +65°C   |       |
| <b>Shock</b>                              | 30g, 11ms, All axes                                     |       |
| <b>Vibration</b>                          | 4.3g, Random vibration, All axes                        |       |
| <b>Humidity</b>                           | 5~95%, no condensation                                  |       |
| <b>Software Support</b>                   |   |       |
| <b>SDK</b>                                | Support   |       |
| <b>PC software</b>                        | Support   |       |

### 3. Installation Guide

#### 3.1 Product and Accessory List

Before installation, please make sure that the equipment in the package is intact and all parts are available.

| Standard packaging list   |  |  |   |
|---|--|--|---|
| <p>Camera</p>                | <p>adapter bracket</p>  | <p>Network cable ×1</p>  | <p>Power cable ×1</p>        |
| <p>Blackbody (1 unit)</p>  | <p>Power adaptor</p>  | <p>Tripod ×2</p>       | <p>Extension shaft ×1</p>  |

#### Optional package (Prepared by users):

- **Computer performance:** Configuration requires i5 and above CPU, 8G and above memory, 64-bit win10 system, supports Gigabit network;
- **Display:** >1920 × 1080 resolution;
- **Power adaptor:** 10 ~ 36VDC, 12VDC is recommended;
- **Tools:** screwdriver, etc.

### 3.2 Mounting and Connecting Thermal Imager

As shown in the figure below, 4 × M2 mounting holes are reserved at the bottom of the product for fixing the camera to the adapter or bracket.

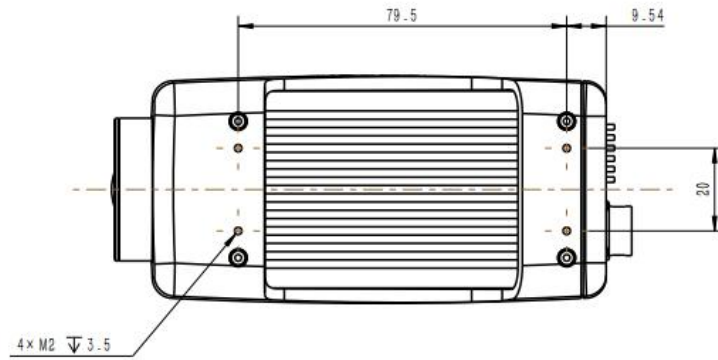


Figure 2 Bottom structure of camera

Take the thermal imager mounted on a tripod as an example, as follows:

- Fix the adapter bracket on the thermal imager with 4 screws(M2\*6), as shown in Figure 3;

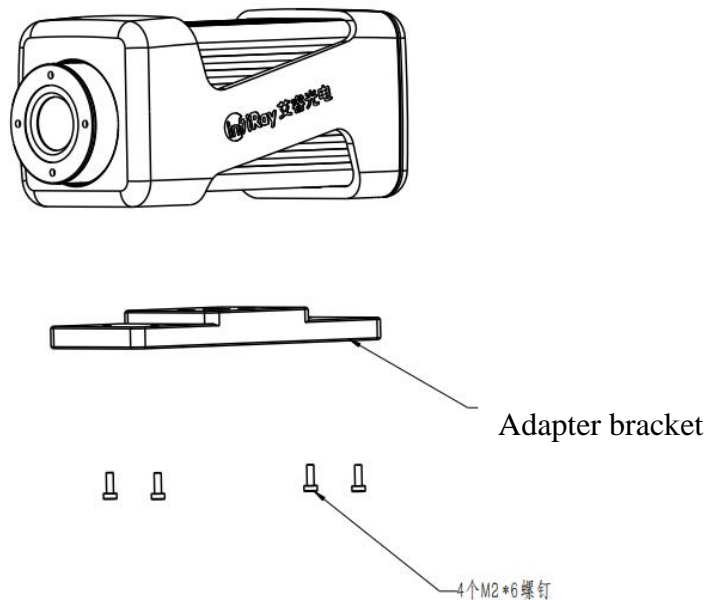


Figure 3 Installation of camera adapter



- After the adapter bracket is installed, the bottom of the thermal imager is shown in Figure 4, and it is mounted on the corresponding supporting equipment through the mounting holes 3 / 8-16UNC or 1 / 4-20UNC at the bottom of the adapter bracket.

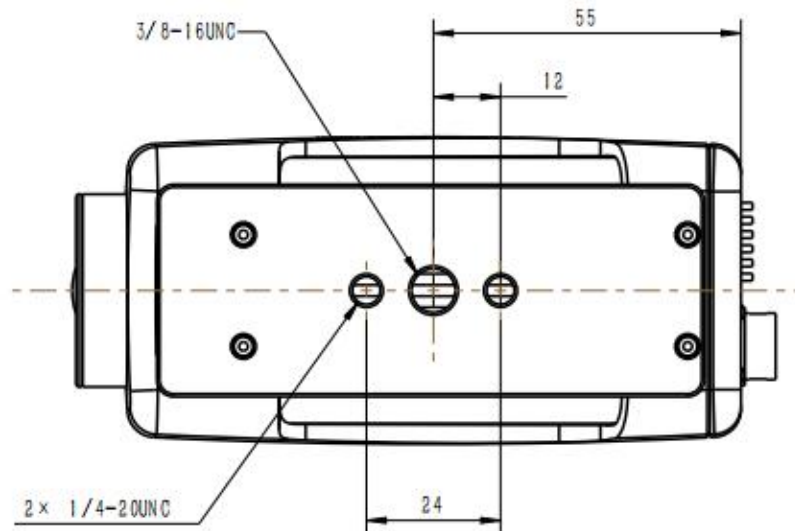


Figure 4 The bottom structure of camera adapter

Take the thermal imager mounting on a tripod as an example, the specific steps are as follows :

Step 1: Remove the fixed adapter on the tripod;

Step 2: Screw the screws on the adapter into the 3 / 8-16UNC mounting holes on camera adapter ;

Step 3: Attach the fixed base to the tripod with the camera.



Step 1



Step 2



Step 3

Figure 5 Thermal camera mounting steps

2) Plug the power cable and network cable into the ports of the camera and the RJ45 interface of the network, the other end of the network cable is connected to the computer network port, and the other end of the power cable is connected to a DC power adaptor (12V is recommended).

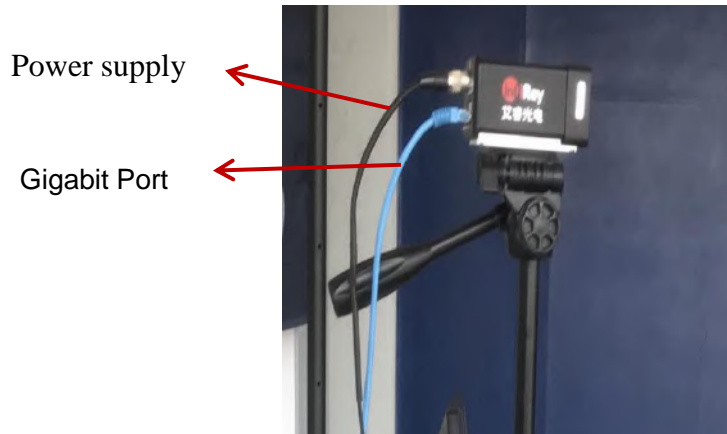


Figure 6 Interface of camera

3) Power up the power adaptor and power up the camera. The sound of the shutter action is heard immediately after power-on, indicating that the camera has been powered on successfully.


4) After the power on 1 minute, run the web client or PC client , connect the device;

5) In order to ensure the accuracy of temperature measurement, the thermal imager needs to be stable for 40 minutes after powering on, and then start the temperature monitoring

**Notes:**

Following picture shows the demension of connector which is used to provide power.

It provides a reference to buy a matched DC power adaptor .

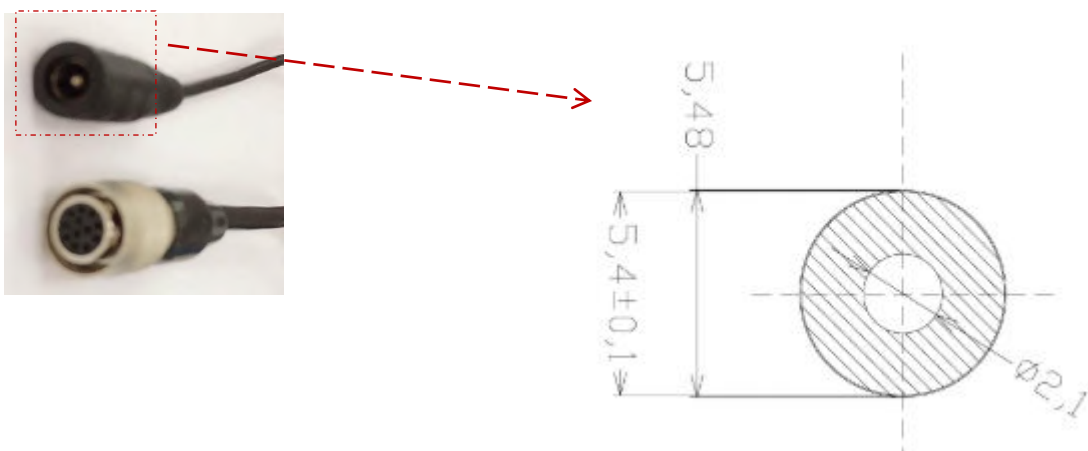

### 3.3 Introduction of Reference Blackbody

#### 3.3.1 Specification

The specification of reference blackbody is shown in table 2.

Table 2 Reference blackbody parameters

| Reference blackbody       |                    |
|---------------------------|--------------------|
| Power supply              | AC180~240V         |
| Max. power consumption    | 260W               |
| Environment temperature   | 0°C~40°C           |
| Accuracy                  | ±0.2°C@35°C        |
| Emissivity                | 0.97±0.02          |
| Settling time             | ≤5min              |
| Temperature stability     | ± (0.1°C~0.3°C) /h |
| Effective radiant surface | 100mm×80mm         |
| Dimension (L×W×H)         | 195mm×150mm×110mm  |

#### 3.3.2 Structure of Reference Blackbody

As shown in figure 7, you can find the drawing of reference blackbody. The top and bottom of the blackbody are reserved with the American standard 1/4-20UNC tripod interface, which is convenient to install for users. Four foot-pads are distributed at the bottom of the four corners, can be used for horizontal placement.

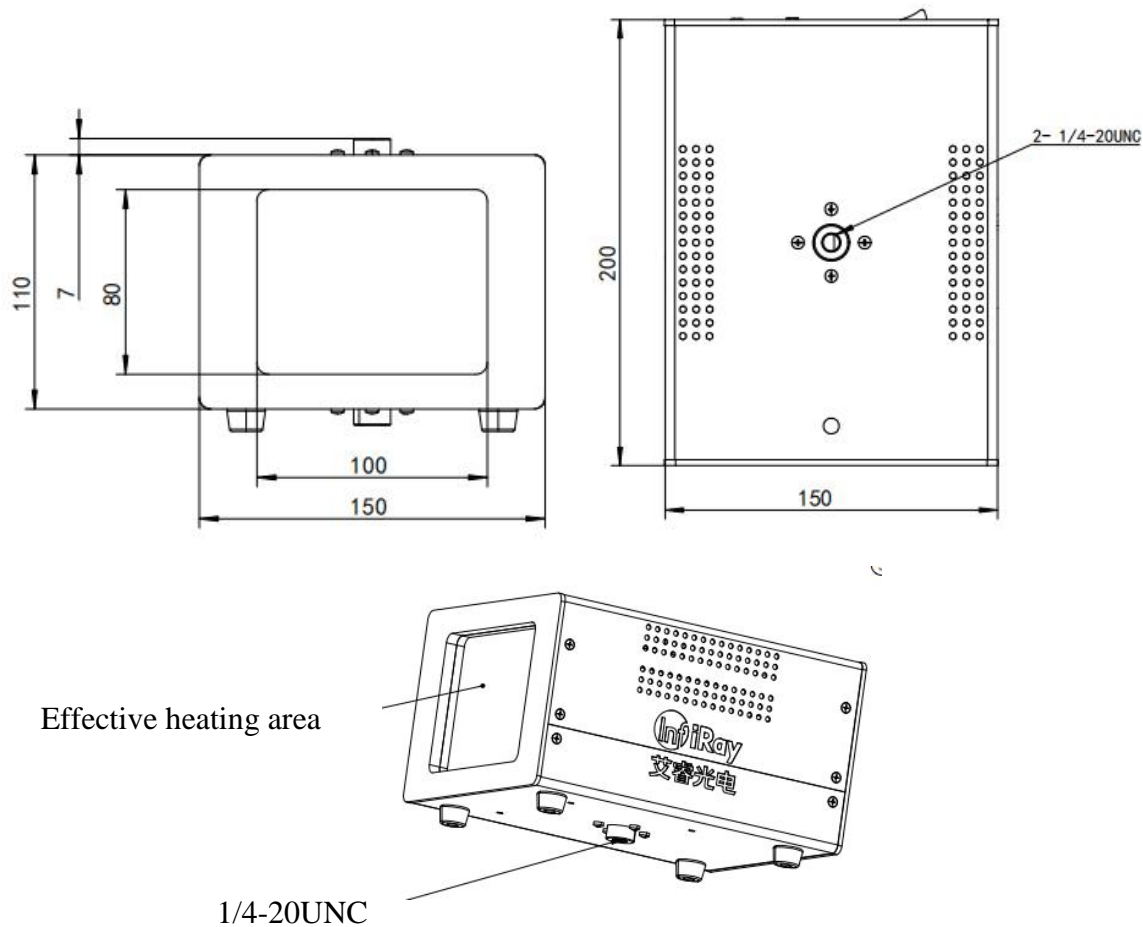


Figure 7 Drawing of blackbody

### 3.3.3 Quick Start and Attentions

Connect power supply(AC180~240V), and turn on the power switcher on the blackbody, the power indicator is on, it display 35°C in default. The blackbody can be used normally when it gets stable after waiting 5 minutes.

**Attention:**

- Please ensure the correct input of power supply;
- Please do not to touch the heating area of blackbody, in case of scalding;
- The blackbody should be used under the environment temperature of 25°C, to ensure the accuracy of equipment;

- Avoid using in the humid environment and beyond the operating temperature range, so as not to affect the accuracy of the equipment.

### 3.4 Site Layout Solution

#### 1) 、 Reference layout

Refer to Figure 8 for the layout plan in a wide space area. The external reference blackbody is placed directly in front of the infrared camera and made it observable in the imaging field of view . Pedestrians need to be guided and make sure all of them can be observed with the thermal imaging camera.

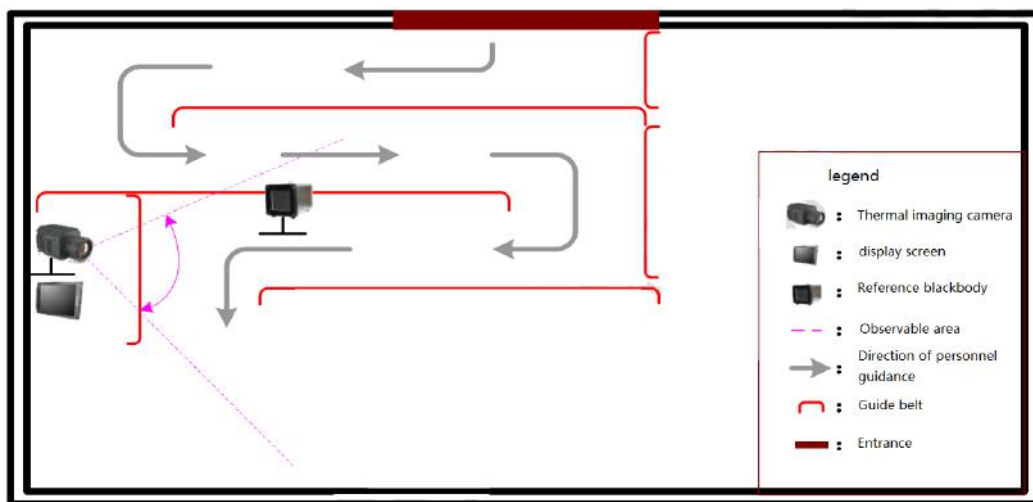


Figure 8 Arrangement of infrared human body thermal imaging camera

- Installation direction and height of camera

The installation distance and height of the camera should be arranged according to the actual situation, the recommended distance is **2m**, and the installation height should ensure that all of the observed pedestrians can appear in the field of view.

- Reference black body layout

As shown in Figure 9, adjust the direction and position of the infrared human body thermal imaging camera and the reference black body so that the reference black body appears in the 2/3 central area of the camera market and shall not be shielded by pedestrians.

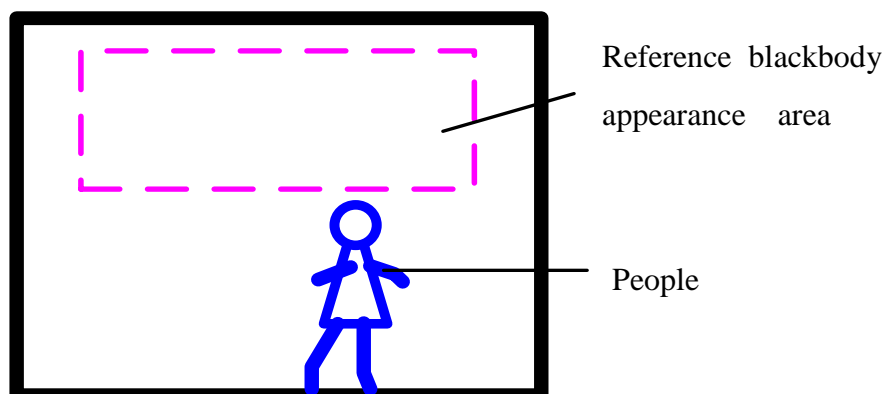


Figure 9 Infrared imaging on PC

- Environmental requirements

The infrared temperature measuring camera is arranged in an environment where the temperature is relatively constant and there is no direct sunlight, to avoid temperature measurement errors caused by the temperature changing environment affecting the human skin surface temperature. According to the above installation requirements, the following installation methods can be used to install equipment in halls and aisles.

### 3.4 Site Installation and Layout Considerations

1) 、 Avoid placing it at the entrance with airflow disturbance

Avoid the temperature rise and fall caused by the air surface disturbance of the skin of the person being tested, which will affect the false alarm or miss alarm of the infrared temperature measuring camera.

2) 、 Select the right monitoring scenario

Avoid using infrared cameras in environments with large areas of high temperature objects (such as electronic billboards, large area light boxes, glass walls that are exposed to the sun, etc.). These hot objects may cause false alarms or missed alarms.

3) 、 The ground of the monitoring site should not be too smooth


Avoid the temperature rise and fall caused by the airflow disturbance on the skin surface of the person being tested, which will affect the false alarm or leakage of the infrared thermometers.

4) 、 Reserve a long enough channel

Prevent the person being tested from entering the room and testing directly in a hot / cold outdoor environment. Since the skin temperature is higher or lower than real skin temperature when they first enter a room with a constant temperature from a hot / cold outdoor environment, at this time, a long enough channel needs to be reserved to allow the subject's body surface temperature to gradually stabilize before testing. To avoid false alarms or missed alarms of the infrared camera.

## 4 Quick Start of PC Client

Professional temperature screening client will be provided , users can easily get the thermal image and body temperature. If you need “InfiRay Intelligent Infrared Temperature Screening System Operation Manual”, please contact us.

| Product model   | tool  | Instructions for using   |
|-----------------|---|--|
| AT3003F/AT6005F |  IRay Intelligent Temperature Screening System | “InfiRay Intelligent Infrared Temperature Screening System Operation Manual” |

please refer to the following method to set environmental parameters on PC client .

- **Atmospheric temperature:** Set according to the ambient temperature of the site.
- **Reflected temperature:** If there is no high temperature object in the use environment, set it to the ambient temperature.
- **Emissivity:** If measuring body temperature, it is recommended to set as 0.98.
- **Atmospheric transmittance:** Use indoors or in places with good visibility. Keep the default parameter 1.
- **Distance:** set according to the use conditions, the default is 2m.
- **Temperature measurement mode:** Standard mode and correction mode are available for selection. The temperature measured in the standard mode is the real temperature. The temperature measured in the correction mode is the internal temperature of the human body calculated from the human skin temperature. Especially, correction1 should be selected when carrying out body temperature screening .

## 5、 Troubleshooting

| Fault description   | Possible Causes   | Method of exclusion   |
|---|---|---|
| Image blur  | Objective focal length does not match   | The thermal imager is aimed at the target position, click the auto focus on the PC client |
|   | Long time without image correction  | On the host computer, click the shutter correction  |
| Boot failure  | The supply voltage exceeds the normal working supply voltage range  | Check if the power supply voltage is between 10 ~ 36V                                     |
|   | Power connector is loose.   | Check whether the power interface is connected reliably.                                  |
| Low measuring accuracy  | The thermal camera stabilization time is too short  | Keep the camera stable for more than 40 minutes   |
|   | Reference blackbody position setting is incorrect   | Reset the position of blackbody   |
|   | Unfocused   | Aim at the target position for auto focus calibration                                     |
| Boot failure with report "without mfc300.dll library"         | No VC environment installed   | Install VC environment library  |
| Click preview and show parameter error                        | Network IP address has not modified to the same IP segment of camera  | Manually configure computer address to 192.168.1.XX network segment                       |
| Network parameters were recognized, but the connection failed | The network IP address is set incorrectly   | Modify the computer's fixed IP address to the correct network segment                     |
| Image stuck   | 1、 Check whether the power and network cable connections are dropped.<br>2、 The power supply is interrupted | After detecting the power, restart the software to connect                                |
| Image video freeze  | Inadequate computer configuration   | Check if the computer support Gigabit Ethernet  |



## 6、 Maintenance

### 6.1 Notice

- Make sure the camera is fixed firmly when the camera is installed on wall, ceiling or pillar.
- The camera shall not fall on the equipment, especially vigorous vibration equipment. Keep the camera away from the presence of magnetic field interference, avoid mounting camera to the surface where occurs vibration or shocks easily. (The camera may be damaged ignoring this item) .
- The camera lens shall not be aimed at the high temperature target such as sun, incandescent lightetc., or it may damage the lens or detector.
- Equipment used indoors should not be installed exposed to rain or very wet areas.
- Do not place the equipment in an environment with corrosive gas, which may cause damage to the equipment.
- Avoid placing equipment in the area where has direct sunlight, poor ventilation , or near heat sources such as heaters (ignoring this can lead to fire hazards).
- Please keep all original packing materials of the equipment properly, so that in case of any problem, use packing materials to pack the equipment and send it to the agent or return it to the manufacturer for treatment. IRay shall not be responsible for any accidental damage in transit caused by materials other than the original packaging.

### 6.2 Cleaning and Maintenance

- Do not use cleaning products that corrode or scratch optical glass components.
- The germanium window surface is coated with anti-reflection coating, dust, grease, fingerprints will produce harmful substances and lead to its performance decline or cause scratches, once found dirt, please follow the following methods.
  - 1) Use a blown balloon or a soft brush to clean the lens surface and avoid dust particles scratching the anti-reflection film on lens surface during the wiping process.
  - 2) Use soft cotton cloth or lens wiping paper to dip in alcohol or lens wiping liquid, gently wipe the lens surface from the middle to the edge, pay attention to do not squeeze the lens surface,

wipe liquid also can not dip in too much, make sure the cloth is wet, but wipe liquid can not be extruded. If still not clean, replace the cloth and repeat operation.

### 6.3 Safety Codes

- Please use standard power supply, to avoid the camera is damaged by over-voltage.
- Short-circuit is prohibited.
- Hitting or beating the product with sharpen object is prohibited ;
- The product shall not be exposed in the high temperature environment which is higher than 80°C.
- Do not put the product into the fire.

## 7. After sales service

The AT300 / AT600 precision body temperature measuring infrared camera is developed by our company and has good after-sales service guarantee such as equipment maintenance and repair. If you have any needs, please contact our company.

## 8. Company information

IRay Technology Co.,Ltd.

Website: [www.IRaytek.com](http://www.IRaytek.com)

Phone: 86-0535-3410623

Fax: 86-0535-3410610

Email: [sales@IRaytek.com](mailto:sales@IRaytek.com)

Address: 11th Guiyang Street, YEDA, Yantai, Shandong Province